

Esquire IVF Workstation is designed with optimized environment for Oocytes & Embryo's

Esquire IVF Workstation can be **customized** as per Embryologist choice

4 Different Models

Basic – 4 feet

Spacious – 5 feet

Dual – 6 feet

ICSI – 6 feet

Monitoring System


Monitor is flushed in the back wall of the working area for an easy view & detailed oocyte identification & grading

Heated Table Working Area

- Aluminium Composed Silicone Heating Technology
- Silicone with stands all environmental conditions across the globe
- A patterned Silicone heating system is composed in aluminium for a better Accuracy & Uniformity

Built in Glass Stage

- Circular or Rectangular ITO based glass stage in the light opening
- 2.5 mm thick glass stage for long durability
- BK 7 grade glass with transparency accuracy 85%

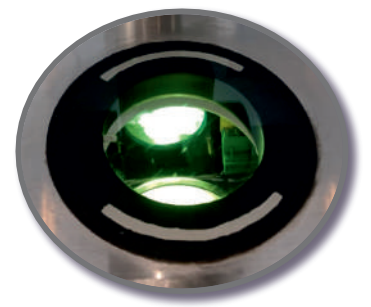

A Product From Your Friend

Esquire IVF Equipments Pvt Ltd

#5A, Sri Mahalakshmi Co-Operative Nagar, Vazhuthalampedu, Kundrathur, Chennai 600 069
www.esquires.in | admin@esquires.in



Microscope Mounting

- Any brand stereo microscope can be mounted on the table
- Microscope is mounted in such a way that the eye piece projects outside the workstation to validate contamination protocol

LED Light Source Transmitted

- LED light source with bi-directional & frosted mirror to focus the light intensity precisely.
- In-built LED yellow/ grey light source for easy embryo access



Miniature Test Tube Warmer

- Built-in miniature test tube warmer movable to the whole working area of the workstation to facilitate the oocyte identification

A Product From Your Friend

Esquire IVF Equipments Pvt Ltd

#5A, Sri Mahalakshmi Co-Operative Nagar, Vazhuthalampedu, Kundrathur, Chennai 600 069
www.esquires.in | admin@esquires.in